

Naída Bimodal Hearing Solution and Roger™: Easy to Connect



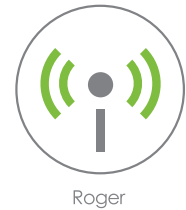
Recipients with the Naída bimodal hearing solution can easily connect to a range of Phonak Roger™ technologies for simultaneous wireless streaming to a Naída CI sound processor and a Phonak Naída™ Link hearing aid.

What is Roger?

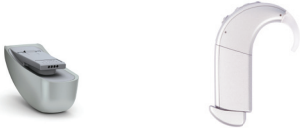
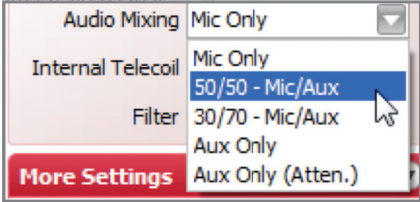

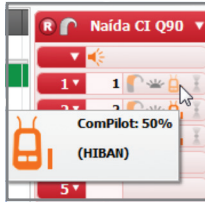

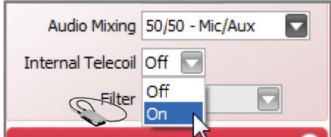
- Roger is the digital wireless standard.
- Roger wirelessly transmits the speaker's voice directly to the listener to improve speech understanding in noise and over distance.
- Roger transmits high-quality audio to Naída CI and Naída Link at the same time.

Compatible Roger Microphone Solutions

- Roger™ Pen
- Roger™ EasyPen
- Roger™ Clip-On Mic
- Roger inspiro™
- Roger Touchscreen Mic
- Roger Pass-around
- Roger Multimedia Hub
- Roger Table Mic



Roger™ and the Naída CI Sound Processor

Options	Roger Receiver for the Naída CI Sound Processor	Fitting Considerations
<p>Roger 17</p>	<p>Roger™ 17 receiver is an integrated receiver that connects seamlessly to the bottom of Naída CI.</p>  <p>Roger 17 attaches to the bottom of the PowerCel™ 170. Roger 17 (02) is compatible with transmitters, including Roger inspiro.</p>	<p>Select an Audio Mixing setting that includes Aux input. (50/50 – Mic/Aux is the default.)</p> 
<p>Roger X</p>	<p>Roger X is a universal miniature receiver that is compatible with Naída CI via the ComPilot accessory.</p>  <p>Roger X receiver attaches to the bottom of the ComPilot.</p>	<p>Select the ComPilot percentage. (50% is the default.)</p> 
<p>Roger MyLink</p>	<p>Roger MyLink is a telecoil neckloop receiver.</p> 	<p>Select an Audio Mixing setting that includes Aux input and enable the Internal Telecoil.</p> 

Roger™ and Naída™ Link Hearing Aid Compatibility

PHONAK

Naída Link UP Hearing Aid



Roger 10 is a design-integrated receiver.



Naída Link RIC Hearing Aid



Roger 15 is a design-integrated receiver.

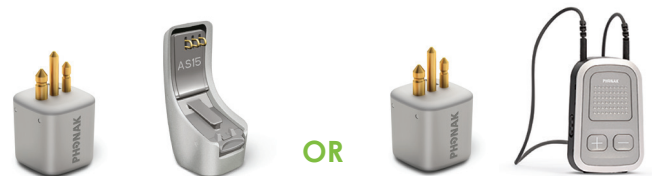


Roger X is a universal miniature receiver that is compatible with Naída Link via the AS10 (audio shoe) or via the ComPilot accessory.



Roger X attaches to the bottom of the audio shoe or to the bottom of the ComPilot.

Roger X is a universal miniature receiver that is compatible with Naída Link via the AS15 (audio shoe) or via the ComPilot accessory.



Roger X attaches to the bottom of the audio shoe or to the bottom of the ComPilot.

Roger MyLink

Roger MyLink is a telecoil neckloop receiver.



Roger and Naída Link Programming Consideration

The Bimodal Fitting Report, generated in SoundWave™, is designed to guide the fitting of Naída Link for bimodal use with Naída CI. The Bimodal Fitting Report provides information on how Naída Link should be programmed for use with Naída CI and Roger.

Roger Support Tools

For additional information on Roger, please visit the Phonak support site for wireless accessories:
<https://www.phonakpro.com/com/en/support/product-support/wireless-accessories.html>